


CLEARPRO Office

Reading 1.25 Shift 1.25

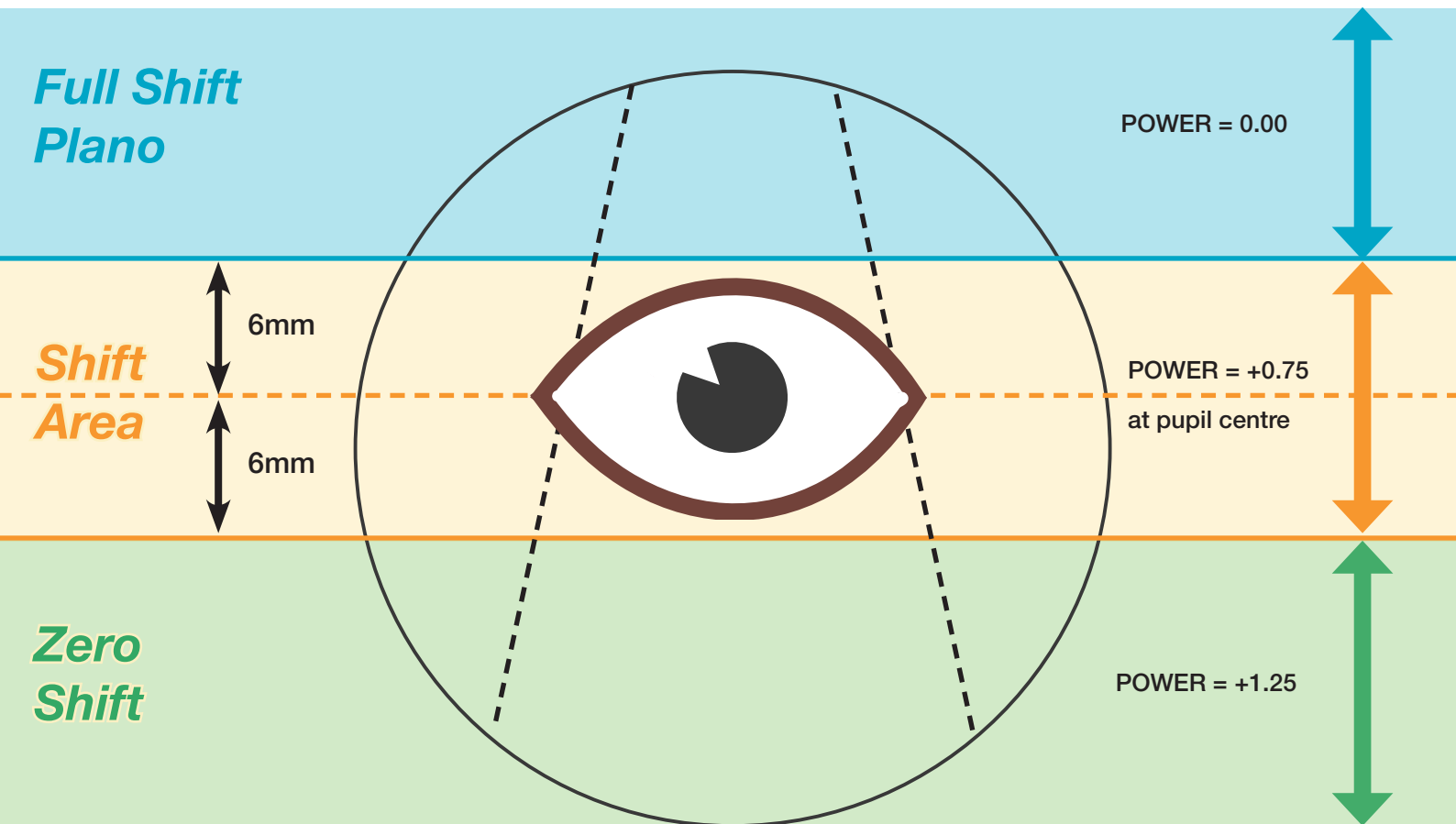


Please Note: Measurements are not to scale

- ORDER AS READING RX
- POWER DECREASES FROM READING AREA
- SHIFT AREA STARTS BELOW LOWER LID
- PUPIL CENTRE = 60% OF SHIFT POWER
- FULL SHIFT STARTS ABOVE UPPER LID

CLEARPRO

Relax

Plano 0.00 Add 1.25

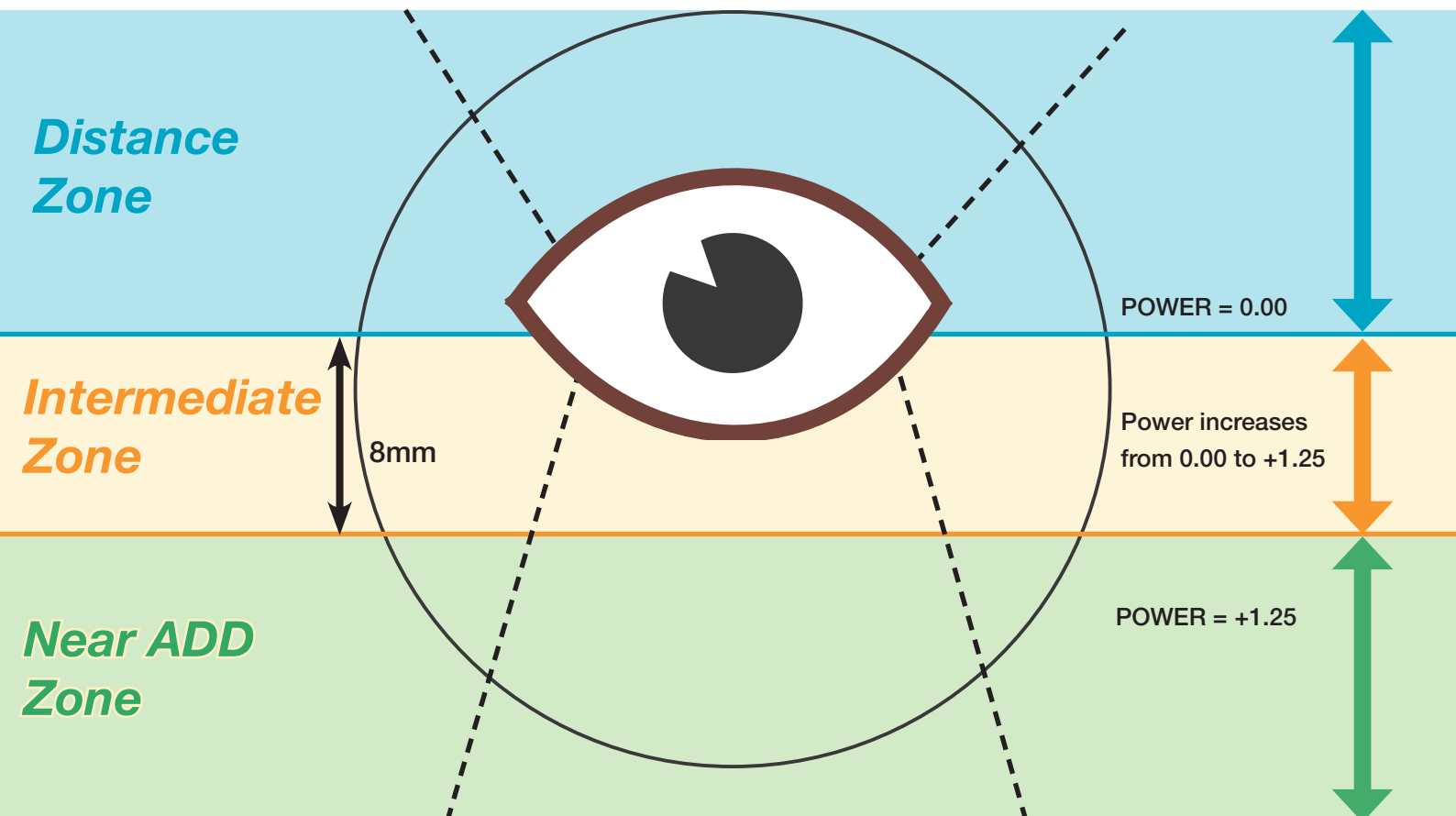


Please Note: Measurements are not to scale

- ORDER AS DISTANCE RX WITH READING ADDITION
- ADD RANGE IS +0.25 TO +2.50
- ORDER THE ADD IN PLUS POWER
- INTERMEDIATE ZONE STARTS 2MM BELOW PUPIL CENTRE
- POWER CONTINUES TO INCREASE BELOW PUPIL HEIGHT

INFINITY 3D MAXIMA

Plano 0.00 Add 1.50

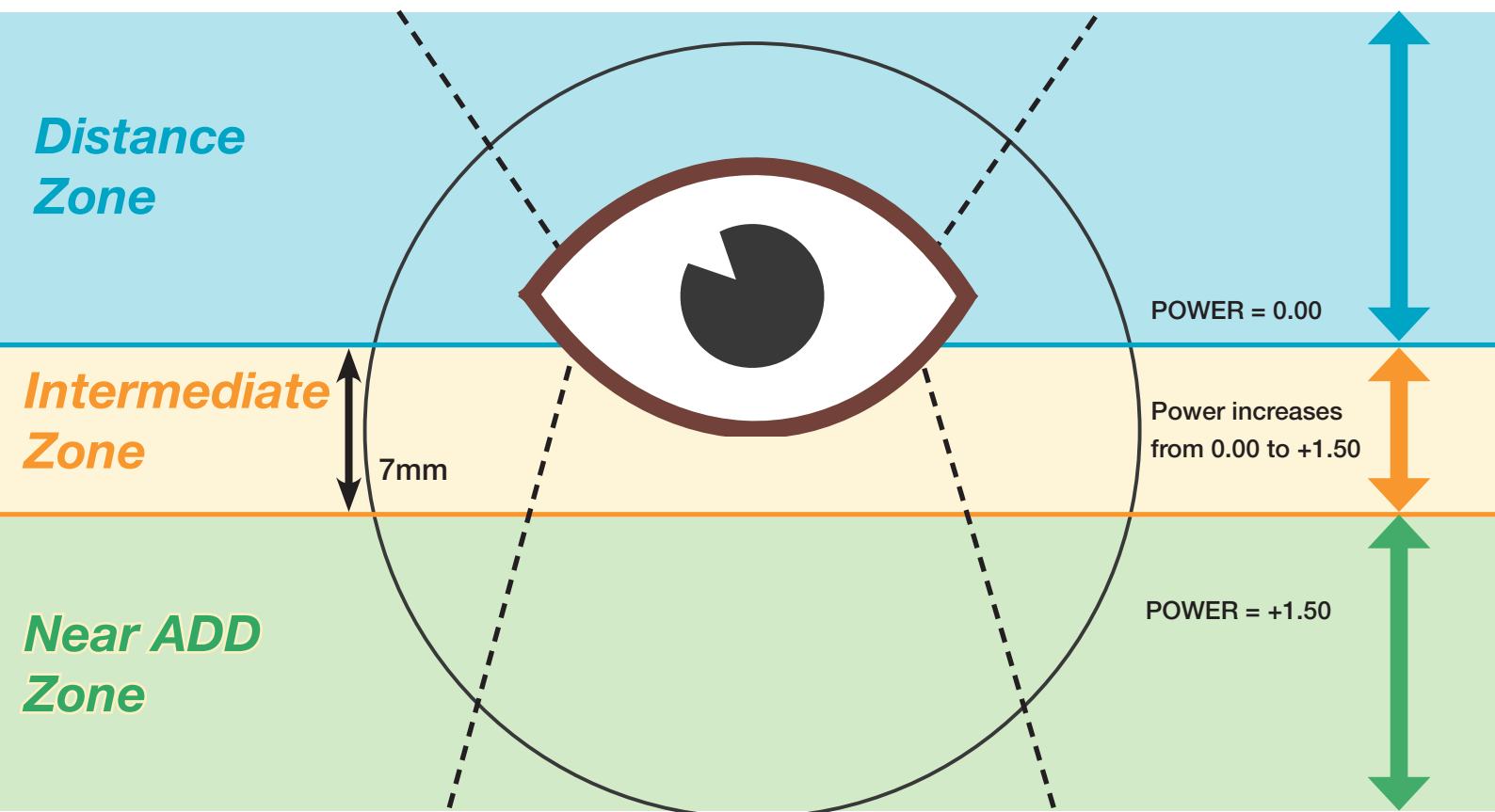


Please Note: Measurements are not to scale

- ORDER AS DISTANCE RX WITH READING ADDITION
- ADD RANGE IS $+0.75$ TO $+3.50$
- ORDER THE ADD IN PLUS POWER
- INTERMEDIATE ZONE STARTS 4mm BELOW PUPIL CENTRE
- POWER CONTINUES TO INCREASE BELOW PUPIL HEIGHT

INFINITY 3D MULTIVIEW

Plano 0.00 Add 1.50

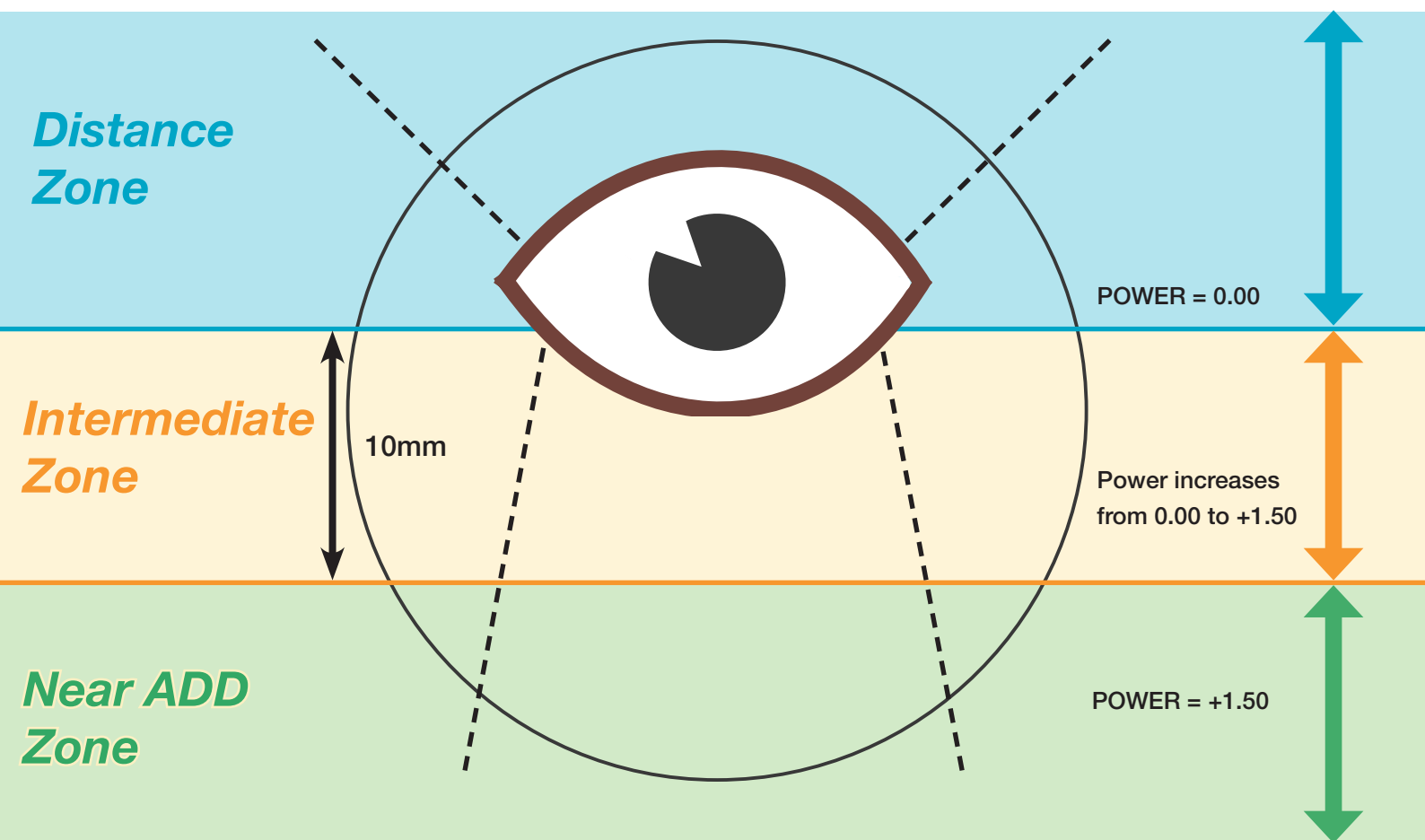


Please Note: Measurements are not to scale

- ORDER AS DISTANCE RX WITH READING ADDITION
- ADD RANGE IS +0.75 TO +3.50
- ORDER THE ADD IN PLUS POWER
- INTERMEDIATE ZONE STARTS 4MM BELOW PUPIL CENTRE
- POWER CONTINUES TO INCREASE BELOW PUPIL HEIGHT

INFINITY 3D OUTDOOR

Plano 0.00 Add 1.50

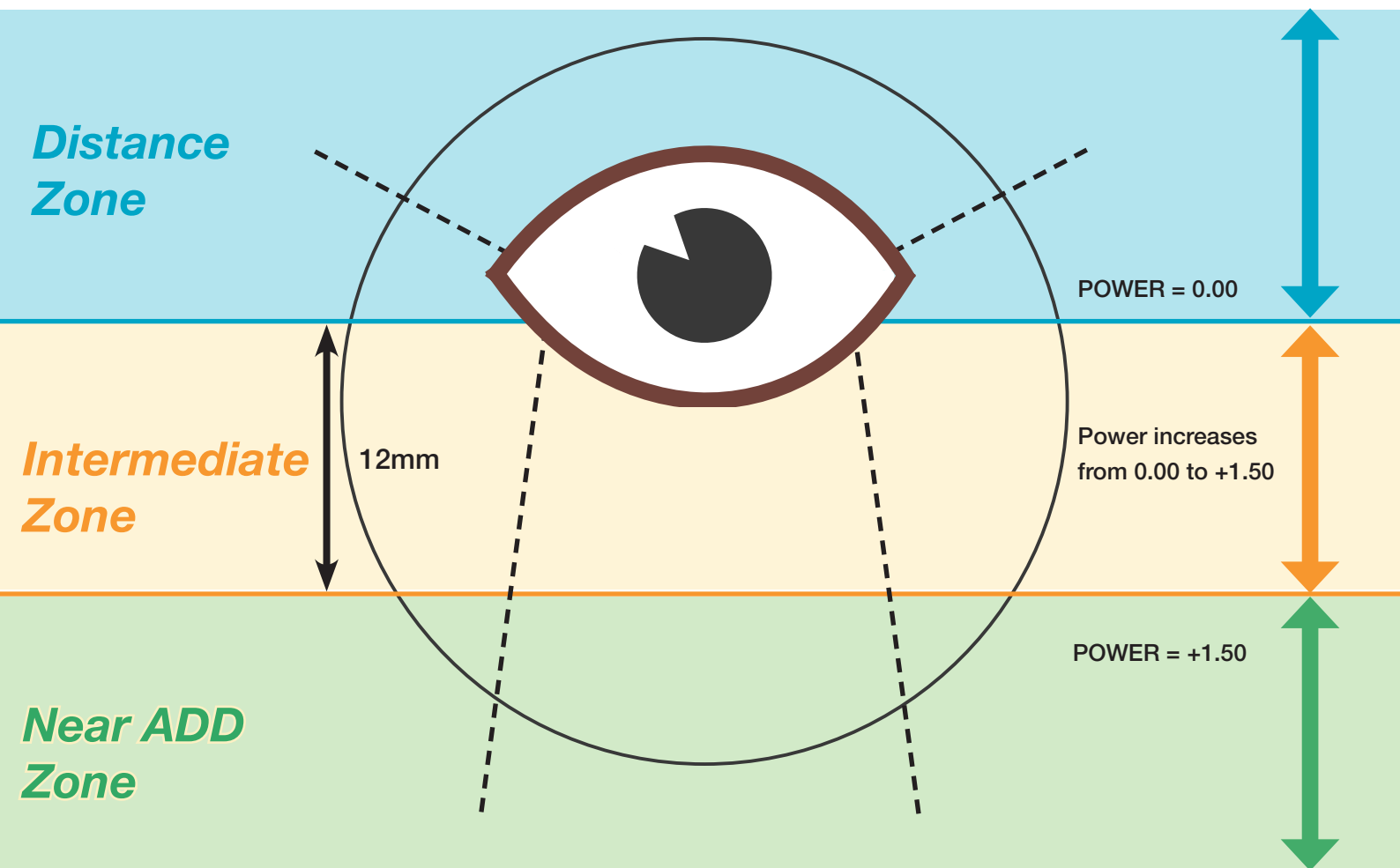


Please Note: Measurements are not to scale

- ORDER AS DISTANCE RX WITH READING ADDITION
- ADD RANGE IS +0.75 TO +3.50
- ORDER THE ADD IN PLUS POWER
- INTERMEDIATE ZONE STARTS 4MM BELOW PUPIL HEIGHT
- POWER CONTINUES TO INCREASE BELOW PUPIL HEIGHT

INFINITY 3D I-ZONE

Plano 0.00 Add 1.50

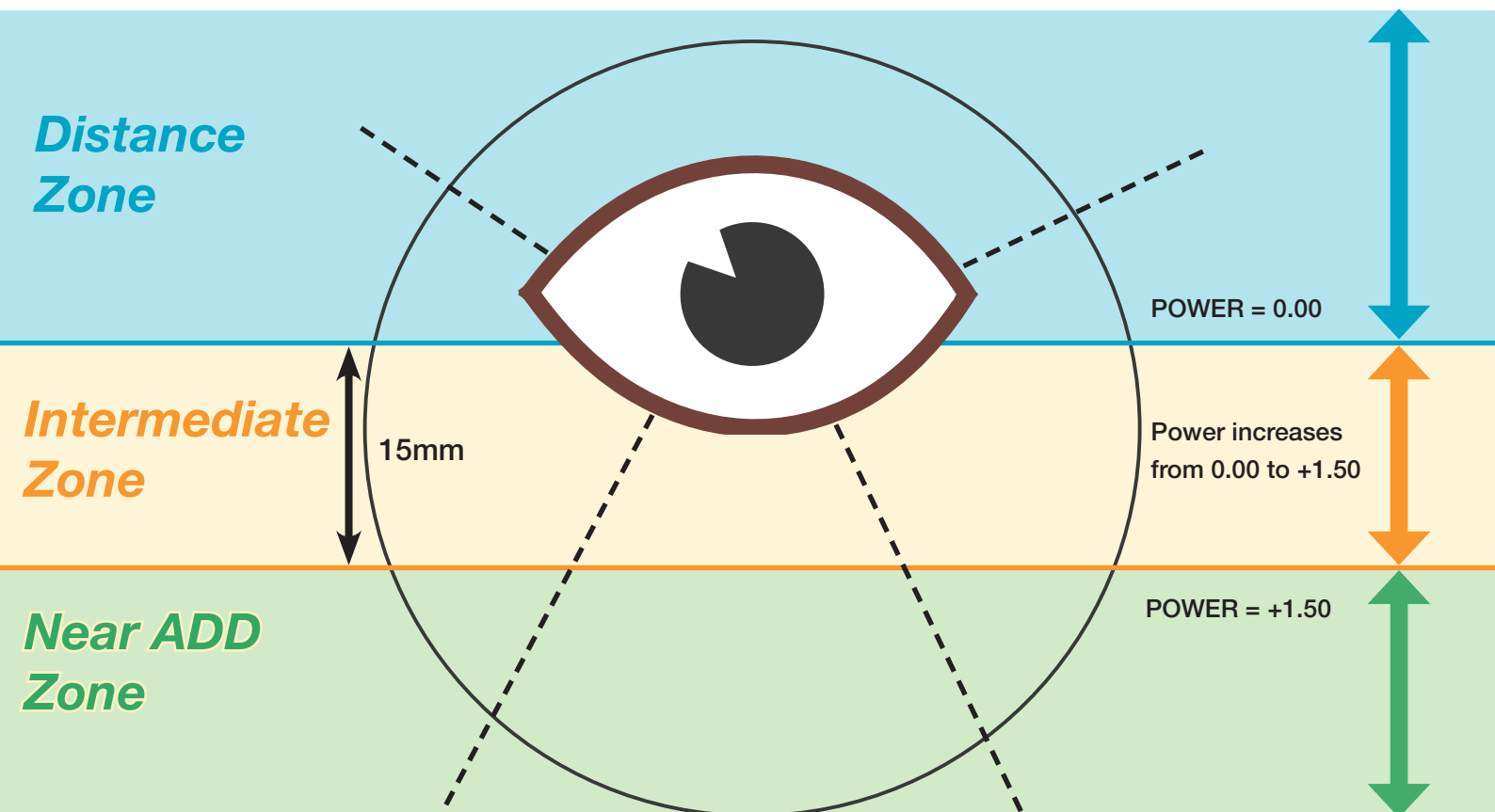

Please Note: Measurements are not to scale

- ORDER AS DISTANCE RX WITH READING ADDITION
- ADD RANGE IS $+0.75$ TO $+3.50$
- ORDER THE ADD IN PLUS POWER
- INTERMEDIATE ZONE STARTS 4mm BELOW PUPIL CENTRE
- POWER CONTINUES TO INCREASE BELOW PUPIL HEIGHT